


STATEMENT OF DEFICIENCIES AND PLAN OF CORRECTIONS	(X1) PROVIDER/SUPPLIER/CLIA IDENTIFICATION NUMBER: 14C0001186	(X2) MULTIPLE CONSTRUCTION A. BUILDING B. WING	(X3) DATE SURVEY COMPLETED 08/02/2022
---	--	--	--

NAME OF PROVIDER OR SUPPLIER WESTERN DIVERSEY SURGICAL CENTER	STREET ADDRESS, CITY, STATE, ZIP CODE 2744 N WESTERN AVE , CHICAGO, Illinois, 60647
--	--

(X4) ID PREFIX TAG	SUMMARY STATEMENT OF DEFICIENCIES (EACH DEFICIENCY MUST BE PRECEDED BY FULL REGULATORY OR LSC IDENTIFYING INFORMATION)	ID PREFIX TAG	PROVIDER'S PLAN OF CORRECTION (EACH CORRECTIVE ACTION SHOULD BE CROSS-REFERENCED TO THE APPROPRIATE DEFICIENCY)	(X5) COMPLETION DATE
Q0000	INITIAL COMMENTS	Q0000		
Q0240	<p>INFECTION CONTROL</p> <p>CFR(s): 416.51</p> <p>The ASC must maintain an infection control program that seeks to minimize infections and communicable diseases.</p> <p>This CONDITION is NOT MET as evidenced by:</p> <p>Based on observation, document review and interview it was determined that the Facility failed to ensure adherence to effective infection control practices related to the high-level disinfection of the equipment used for vaginal ultrasounds and the ongoing implementation of the infection control surveillance of post-surgical infections. This potentially affects all patients that present to the Facility that require surgical services. As a result, the Condition for Coverage, 416.51 Infection Control was not in compliance.</p> <p>Findings include:</p> <p>The Facility failed to ensure that the infection control program included post-procedural infection surveillance and detection through ongoing data collection and analysis. See deficiency cited at Q-0242. The Facility failed to ensure ongoing collection of data of post procedure infections was included as an integral part of their quality assessment and performance improvement program to ensure ongoing surveillance of post procedure infections. See deficiency cited at Q-0244. The Facility failed to ensure that an infection control plan/program was fully implemented by failing to ensure to follow the manufacturer's</p>	Q0240		

Any deficiency statement ending with an asterisk (*) denotes a deficiency which the institution may be excused from correcting providing it is determined that other safeguards provide sufficient protection to the patients. (See reverse for further instructions.) Except for nursing homes, the findings stated above are disclosable 90 days following the date of survey whether or not a plan of correction is provided. For nursing homes, the above findings and plans of correction are disclosable 14 days following the date these documents are made available to the facility. If deficiencies are cited, an approved plan of correction is requisite to continued program participation.

LABORATORY DIRECTOR'S OR PROVIDER/SUPPLIER REPRESENTATIVE'S SIGNATURE 	TITLE Facility Administrator	(X6) DATE 8/16/2022
--	---------------------------------	------------------------

STATEMENT OF DEFICIENCIES AND PLAN OF CORRECTIONS	(X1) PROVIDER/SUPPLIER/CLIA IDENTIFICATION NUMBER: 14C0001185	(X2) MULTIPLE CONSTRUCTION A. BUILDING B. WING	(X3) DATE SURVEY COMPLETE 08/02/2022
---	--	--	---

NAME OF PROVIDER OR SUPPLIER WESTERN DIVERSEY SURGICAL CENTER	STREET ADDRESS, CITY, STATE, ZIP CODE 2744 N WESTERN AVE , CHICAGO, Illinois, 60647
--	--

(X4) ID PREFIX TAG	SUMMARY STATEMENT OF DEFICIENCIES (EACH DEFICIENCY MUST BE PRECEDED BY FULL REGULATORY OR LSC IDENTIFYING INFORMATION)	ID PREFIX TAG	PROVIDER'S PLAN OF CORRECTION (EACH CORRECTIVE ACTION SHOULD BE CROSS-REFERENCED TO THE APPROPRIATE DEFICIENCY)	(X5) COMPLETE DATE
Q0240	Continued from page 1 instruction for use to perform high-level disinfection, on the ultrasound machine used to conduct vaginal ultrasounds. See deficiency cited at Q-0245.	Q0240		
Q0242	<p>INFECTION CONTROL PROGRAM</p> <p>CFR(s): 416.51(b)</p> <p>The ASC must maintain an ongoing program designed to prevent, control, and investigate infections and communicable diseases. In addition, the infection control and prevent program must include documentation that the ASC has considered, selected, and implemented nationally recognized infection control guidelines.</p> <p>This STANDARD is NOT MET as evidenced by:</p> <p>Based on document review and interview, it was determined that the Facility's failed to ensure that the infection control program included post-procedural infection surveillance and detection through ongoing data collection and analysis.</p> <p>Findings include:</p> <p>1. On 08/01/2022 the Facility's policy titled, "Infection Control Plan" dated 03/01/2008 was reviewed and included, "...The infection control plan will be monitored and evaluated by the Performance Committee ...1. Infection control data will be collected, analyzed and trended. Information obtained will be given to the Office Manager or designee, and used to improve patient care, as well as improve practice's performance in the implementation of its infection/exposure control plan ..." The Facility's was requested to provide the infection control surveillance of data tracking and trending for a sample of selected patients. The Facility was unable to provide documents regarding infection control surveillance. 2. On 08/02/2022 at 11:30 AM, the Nurse Manager (E #2) an interview was conducted. E #2 stated that there is an infection control plan, but there is no data being collected if any patients developed infections after the procedure. E #2 stated that the Nurse Practitioner (E #3) makes the phone calls and is not sure if the nurse practitioner obtains information on patient for data collection of post-procedure infections. 3. On 08/02/2022 at 9:20 AM, the Nurse Practitioner (E #3) was interviewed. E #3 stated that there is</p>	Q0242		

STATEMENT OF DEFICIENCIES AND PLAN OF CORRECTIONS	(X1) PROVIDER/SUPPLIER/CLIA IDENTIFICATION NUMBER: 14C0001185	(X2) MULTIPLE CONSTRUCTION A. BUILDING B. WING	(X3) DATE SURVEY COMPLETED 08/02/2022
---	--	--	--

NAME OF PROVIDER OR SUPPLIER WESTERN DIVERSEY SURGICAL CENTER	STREET ADDRESS, CITY, STATE, ZIP CODE 2744 N WESTERN AVE , CHICAGO, Illinois, 60647
--	--

(X4) ID PREFIX TAG	SUMMARY STATEMENT OF DEFICIENCIES (EACH DEFICIENCY MUST BE PRECEDED BY FULL REGULATORY OR LSC IDENTIFYING INFORMATION)	ID PREFIX TAG	PROVIDER'S PLAN OF CORRECTION (EACH CORRECTIVE ACTION SHOULD BE CROSS-REFERENCED TO THE APPROPRIATE DEFICIENCY)	(X5) COMPLETION DATE
Q0242	Continued from page 2 five percentage possibility of patient developing infection after a surgical procedure. E #3 stated that there no data collected from the patients of post procedural infections. E #3 stated "We do not do any type of surveillance after the procedure." 4. On 08/02/2022 at 10:00 AM, the Medical Director (MD #1) was interviewed. MD #1 stated that he was not aware that the Facility was not conducting surveillance of post procedural infections. MD #1 stated that there is risk for infection after any surgical procedure. MD #1 stated that they plan to hire a consultant to oversee the infection control practices at the facility.	Q0242		
Q0244	INFECTION CONTROL PROGRAM - QAPI CFR(s): 416.51(b)(2) [The program is -] An integral part of the ASC's quality assessment and performance improvement program This STANDARD is NOT MET as evidenced by: Based on document review and interview, it was determined that the Facility failed to ensure ongoing collection of data of post procedure infections was included as an integral part of their quality assessment and performance improvement program to ensure ongoing surveillance of post procedure infections. Findings include: 1. On 08/01/2022, the Consulting Committee Meeting Minutes for 01/2021 – 07/2022 was reviewed. The meetings were held quarterly, and sign-in sheet included the key members of the facility attending the committee meetings. The committee meeting minutes included topics of discussion such as: credentialing, approval of policies and procedures, tissue review report, the infection control report on hand hygiene, infection control prioritized risk for infection control -COVID pandemic, environment of care, census report and employee related agenda. The Consulting Committee Meeting minutes did not include data of the ongoing surveillance of post-procedure infections. 2. On 08/02/2022 at 11:30 AM, the Nurse Manager (E #2) was interviewed. E #2 stated that there is an infection control plan, but there is no data being collected if any patients developed infections after the procedure. E #2 stated that there is no	Q0244		

STATEMENT OF DEFICIENCIES AND PLAN OF CORRECTIONS	(X1) PROVIDER/SUPPLIER/CLIA IDENTIFICATION NUMBER: 14C0001185	(X2) MULTIPLE CONSTRUCTION A. BUILDING B. WING	(X3) DATE SURVEY COMPLETE 08/02/2022
---	--	--	---

NAME OF PROVIDER OR SUPPLIER WESTERN DIVERSEY SURGICAL CENTER	STREET ADDRESS, CITY, STATE, ZIP CODE 2744 N WESTERN AVE , CHICAGO, Illinois, 60647
--	--

(X4) ID PREFIX TAG	SUMMARY STATEMENT OF DEFICIENCIES (EACH DEFICIENCY MUST BE PRECEDED BY FULL REGULATORY OR LSC IDENTIFYING INFORMATION)	ID PREFIX TAG	PROVIDER'S PLAN OF CORRECTION (EACH CORRECTIVE ACTION SHOULD BE CROSS-REFERENCED TO THE APPROPRIATE DEFICIENCY)	(X5) COMPLETE DATE
Q0244	Continued from page 3 data that is presented regarding patients developing infections after the procedure at the consulting committee meetings. 3. On 8/02/2022 at 10:00 AM, the findings were discussed with the Medical Director (MD #1). MD #1 stated that he is not aware of any data regarding patient's developing any type of infections being presented at the consulting committee meetings. MD #1 stated that they plan to hire a consultant to oversee the infection control practices at the	Q0244		
Q0245	<p>INFECTION CONTROL PROGRAM</p> <p>CFR(s): 416.61(b)(3)</p> <p>The program is -</p> <p>Responsible for providing a plan of action for preventing, identifying, and managing infections and communicable diseases and for immediately implementing corrective and preventive measures that result in improvement.</p> <p>This STANDARD is NOT MET as evidenced by:</p> <p>Based on document review and interview it was determined that the Facility failed to ensure that an infection control plan/program was fully implemented by failing to ensure to follow the manufacturer's instruction for use to perform high-level disinfection, on the ultrasound machine used to conduct vaginal ultrasounds. This potentially placed all patients that required a vaginal ultrasound at risk for infection and communicable diseases.</p> <p>Findings include:</p> <p>1. On 08/01/2022 at approximately 10:00 AM, the Facility's policy titled, "Cleaning, Sterilization, and Disinfection Guidelines" dated 02/20/2018 was reviewed and included, "...A.3: The sterile processing technician will review the instruments manufacturer's IFU (instructions for use) to determine the requirement for replicating the validated cleaning and processing methods ...The manufacturer's IFU will be reviewed periodically and processing practices will comply with the most current IFU ...Surgical instruments and equipment should be cleaned and decontaminated according to the manufacturer's validated, written IFU ...methods used for disinfection will be based on ...semi-critical ...risk of infection to the patient ...semi-critical: items that come in contact with</p>	Q0245		

American Health Care Centers

DBA: Western-Diversey Surgical Center

2744 N. Western Ave., Chicago, IL 60647

Dept of Health & Human Services
CMS – Chicago, Survey and Operation Group
233 North Michigan Ave., Suite 600
Chicago, IL 60601-5519

Attn: Wilda Melendez/Annette Hodge

August 17, 2022

We at Western-Diversey Surgical Center are planning to resume procedures on August 24, 2022 (Wednesday).

Sincerely;

A large black rectangular redaction box covering the signature area.

J. Swanson

Administrator

IDPH RESPONSE FORMAT

Prefix Tag Paragraph		Title	Plan of Corrective Action	Date
Q0240	4	Surveillance of post-procedure infections	<p>Because of the sensitivity of the procedures being performed, the consulting committee has decided on writing a letter in a pre-formatted form to survey all patients that had procedures performed to follow-up for post-operative infections. This data will be collected by the infection control coordinator, analyzed and trended weekly for the next 4 weeks and quarterly thereafter and will be reported to the infection control meeting which is part of the quarterly consulting meeting. The collection of data will be started on Sep. 2022. Please see attached format of follow-up form.</p>	08/16/22
Q-0244		HLD	<p>All procedures at the facility was temporarily postponed until deficiency has been resolved.</p> <p>The consulting committee has approved the purchase of Trophone EPR for High Level Disinfection of ultrasound probes.</p> <p>Staff took online course with certificate in preparation for the Trophon training.</p> <p>Staff training on the use of Trophon EPR was conducted by Nanosonics field trainer for the use of Trophon EPR</p> <p>3 of the staff was also trained to be "Super Trainer" for the annual competency of the staff.</p> <p>The staff will be monitored for comptency on a monthly basis for the next 6 months and annully therefater.</p>	07/27/22
				08/03/22
				08/13/22
				08/16/22
				08/16/22

Q-0242

Infection Control Plan

<p>The Infection Control Plan was revisited by the infection control committee and a plan to prioritize risk for infection control was done. The prioritized risk was post-operative infection follow-up of patients. For the purpose of Infection Control follow-up of patients: a letter was created to be sent to the provider to conduct surveillance of post-operative infection of the patient when they follow-up with their provider, the data will be collected monthly, analyzed, trended and will be reported to the Infection Control Committee which meets quarterly with the consulting committee.</p>	<p>08/16/22</p>
--	-----------------

Q-0245

Infection Control Program
"Cleaning and Sterilization
of Instruments"

<p>All procedures at the facility was temporarily postponed until deficiency has been resolved.</p>	<p>07/27/22</p>
<p>The consulting committee has approved the purchase of Trophone EPR for High Level Disinfection of ultrasound probes.</p>	<p>08/03/22</p>
<p>Staff took online course with certificate in preparation for the Trophon training.</p>	<p>08/13/22</p>
<p>Staff training on the use of Trophon EPR was conducted by Nanosonics field trainer for the use of Trophon EPR</p>	<p>08/16/22</p>
<p>3 of the staff was also trained to be "Super Trainer" for the annual competency of the staff.</p>	<p>08/16/22</p>
<p>Policy and Procedure for High Level Disinfection of Ultrasound Probes was created and approved by the consulting committee for implementation Addendum: The log book for the Trophon' EPR HLD will be monitored daily for the next 3 months, and monthly thereafter by the Infection Control Coordinator.</p>	<p>08/16/22</p>

Q-0242


Follow-up for Infection

Q-0244

<p>A post op infection follow-up for elective abortion will be conducted and data will be collected for patients that follow-up or were able to call. Data will be collected by the nurse practioner on a monthly basis, analyzed, trended and reported to the infection control committee.</p>	<p>8/16/2022</p>
---	------------------

Q-0245

<p>The Infection Control Plan was revisited by the Infection Control Committee and prioritized risks was identified, for which a planned data collection was presented. The surveillance indicators planned were: Hand Hygiene, Personal Protective Equipment (PPE), Injection Safety and Medication Handling, Equipment Reprocessing, Environmental Cleaning and Handling og BG minitoring equipt. The above indicators data collection will be conducted daily during procedure days, collected monthly, analyzed and trended and will be reported to the quarterly infection control committee meeting. These activities will be conducted for the next 3 months on daily basis, weekly thereafter and monthly thereafter.</p> <p>Please see attached sample form.</p>	8/17/2022
<p>Please see attached supprting document that the 2 probes with corresponding SN are subjected successfully to Trophon HLD.</p>	8/17/2022
<p>In-service was conducted with the staff on cleaning and disinfection of surgical equip. with Spaulding Classification, in conjunction Trophon EPR was instituted for use after extensive training and certification of staff the manufacturers instructions will be followed on the Trophone EPR for HLD. Oversight will be conducted by the Infection control coordinator on a daily basis for the next 3 months and monthly thereafter</p>	

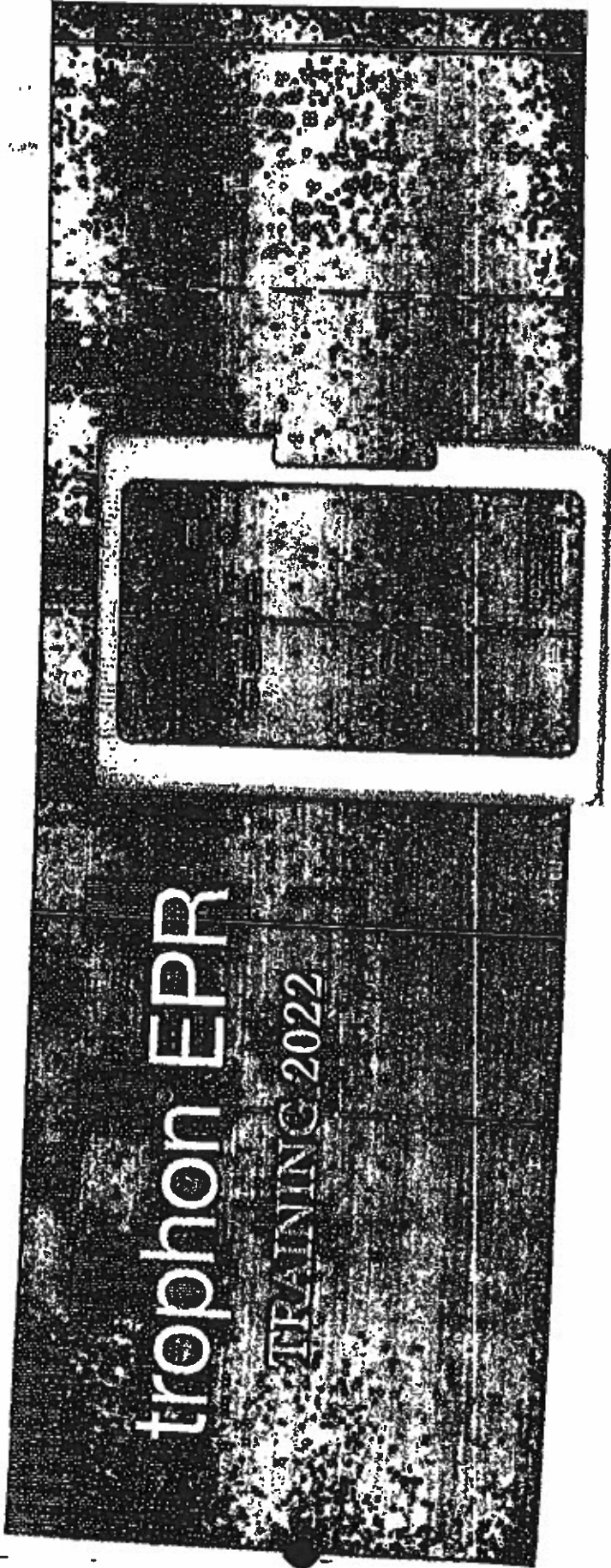


This certifies that

Leslie Aguirre

has successfully completed the trophon EPR product training approved by Nanosonics Limited.
This certificate expires 12 months after the date of issue.

Date of issue: 10 Aug 2022

nanosonics
Infection Prevention For Life

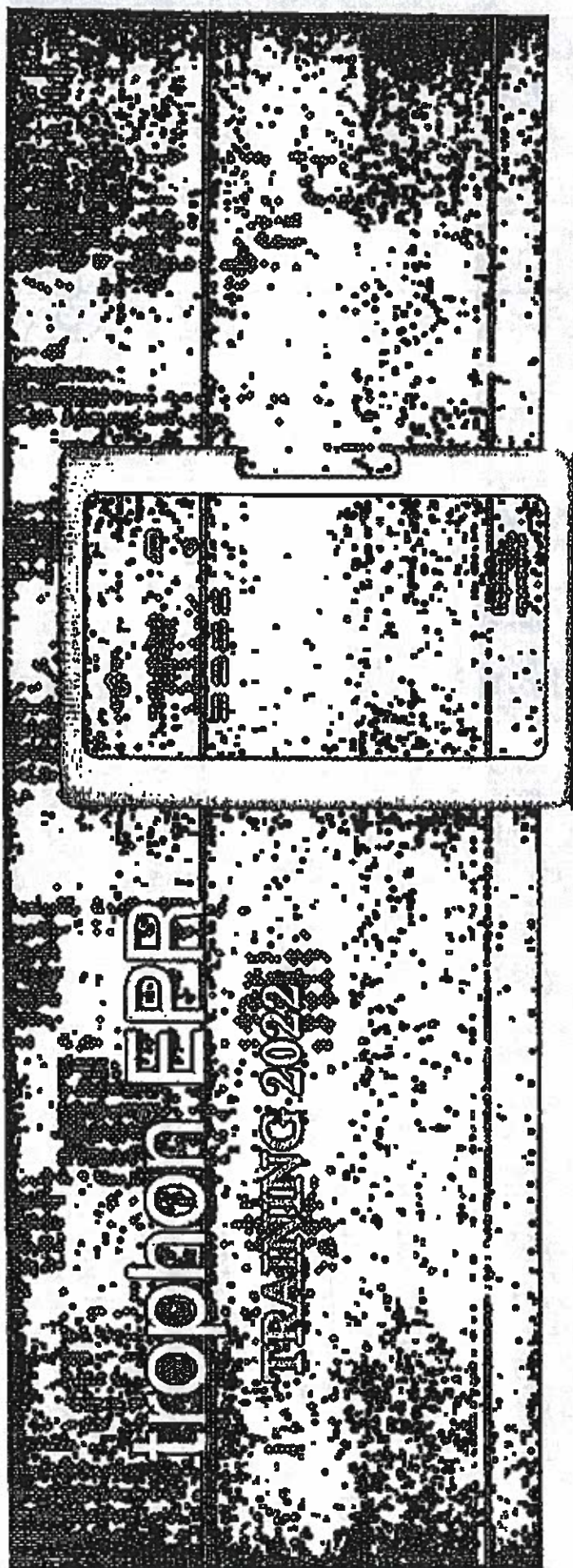

This certifies that


Alexa Tulea

has successfully completed the trophon EPR product training approved by Nanosonics Limited.

This certificate expires 12 months after the date of issue.

Date of issue: 10 Aug 2022



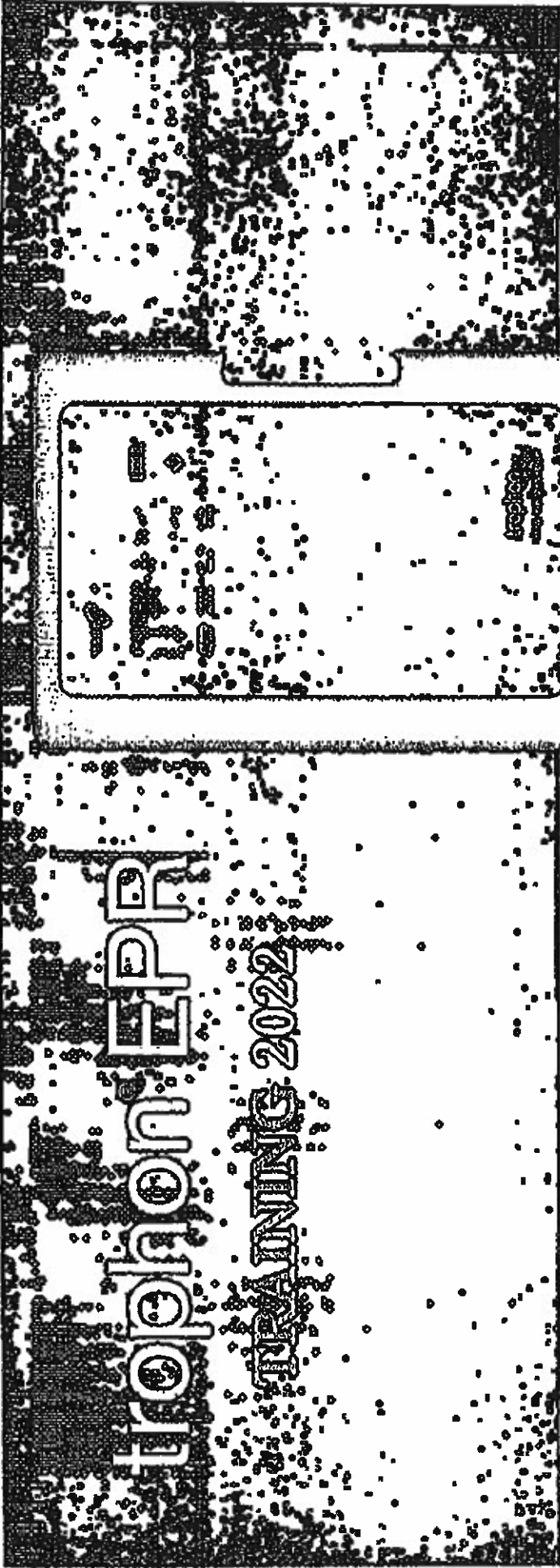

This certifies that


Liauw Devi Gunawan

has successfully completed the trophon EPR product training approved by Nanosonics Limited.

This certificate expires 12 months after the date of issue.

Date of issue: 10 Aug 2022

nanosonics
 Infection Prevention. For Life.

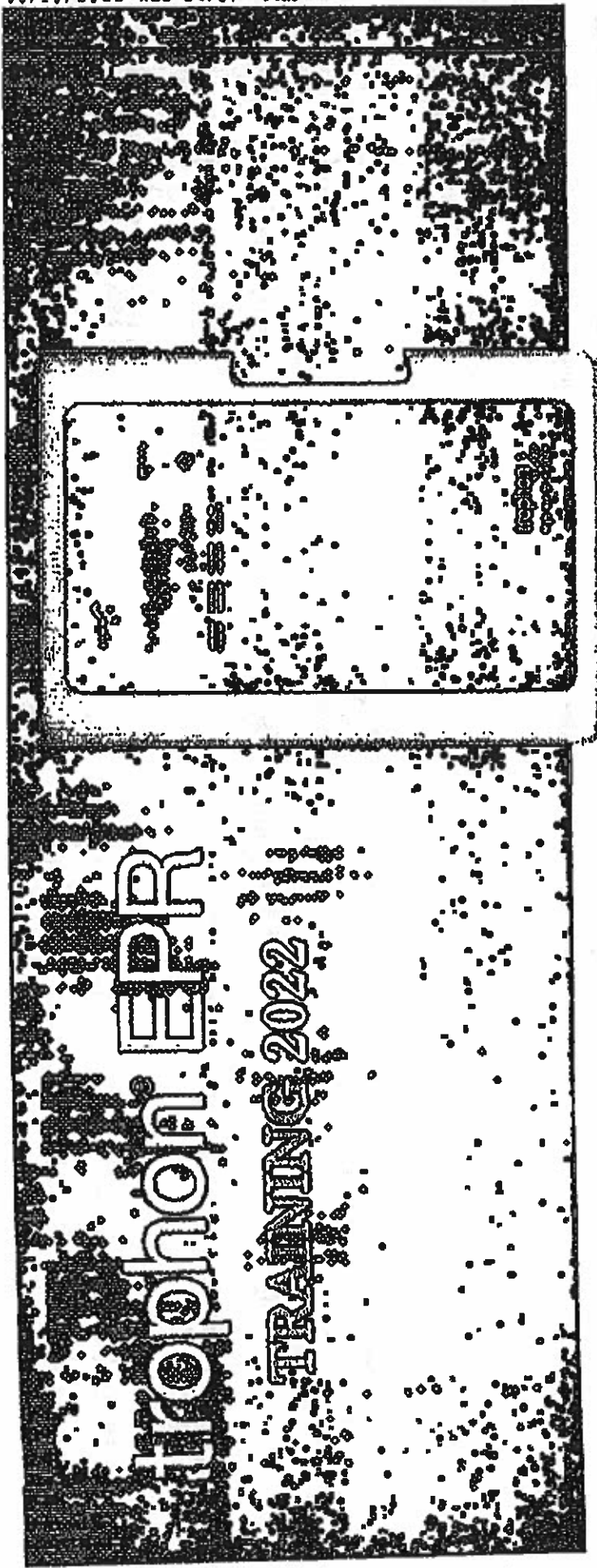

This certifies that

Marie Frukacz

has successfully completed the trophon EPR product training approved by Nanosonics Limited.

This certificate expires 12 months after the date of issue.

Date of issue: 10 Aug 2022

nanosonics
 Infection Prevention For Life.



This certifies that


Leslie Aguirre

has successfully completed the trophon EPR product training approved by Nanosonics Limited.

This certificate expires 12 months after the date of issue.

Date of issue: 10 Aug 2022



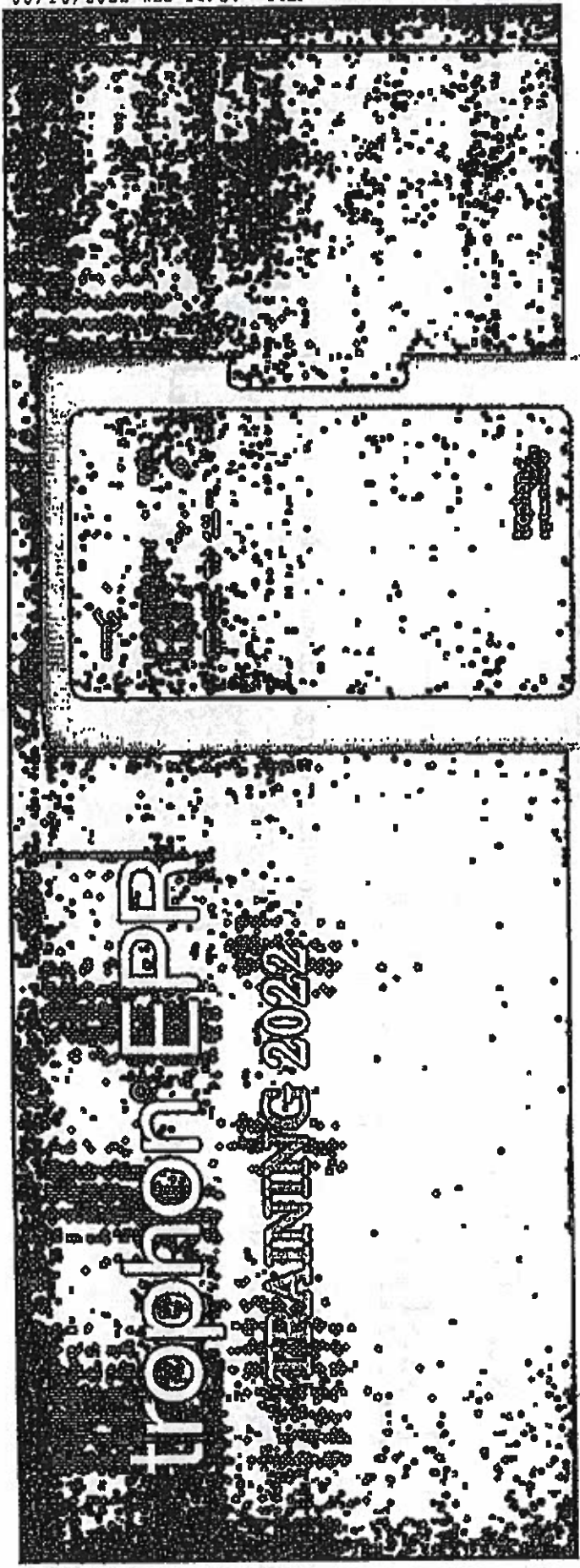


This certifies that

NATALIYA KUKURUZA

has successfully completed the trophon EPR product training approved by Nanosonics Limited.

This certificate expires 12 months after the date of issue.

Date of issue: 10 Aug 2022

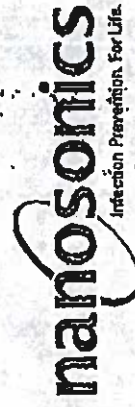

Appendix 6: TROPHON® COMPETENCY SHEET

NAME Marie Fukacz DATE 8/18/22

1. Put gloves on. Note: wear gloves at all times when required as described in the Trophon EPR user manual.

MET NOT MET

2. Pre-clean the probe before the high-level disinfection (HLD) cycle, following the probe manufacturer's instructions for use (IFUs).

MET NOT MET

3. Ensure the probe is clean and free of all visible debris, bioburden, gel, or other soil.
- a. Use Sani Cloths to clean
 - b. Wipe the transducer cord and all surfaces of the transducer until it is visually clean. Use friction and work from cleanest to dirtiest areas.
 - c. Dry the transducer with a soft, dry cloth.
 - d. Visually inspect the transducer to ensure it is both clean and dry prior to HLD.

MET NOT MET

4. Load the clean, dry probe into the Trophon disinfection chamber ensuring:
- a. The probe is secured high in the chamber with tip of probe above embossed line.
 - b. Probe does not contact the chamber wall at any point.

MET NOT MET

5. Place a new red Trophon chemical indicator (CI) into the indicator holder with the red side facing up.
- a. Note: a new CI is to be used for every cycle.

MET NOT MET

6. Close the chamber door and confirm whether the probe is both clean and dry.
- a. If yes, press Start.
 - b. If no, follow the LCD screen prompts.

MET NOT MET

7. At the end of the 7 minute HLD cycle, Trophon's LCD screen states: "CYCLE COMPLETE REMOVE AND WIPE PROBE."

MET NOT MET

8. Open the chamber door.

MET NOT MET

Cleaning, Disinfection, and Sterilization of Patient-Care Items

9. Remove CI, check CI color against the color chart on the CI carton and discard.

a. Note: BOTH the CI and LCD screen must indicate a successful cycle for the probe to be ready for use. If either the CI or Trophon LCD screen indicates a fail, repeat the cycle.

MET NOT MET

10. Remove and wipe the probe using a clean, dry single use cloth.

MET NOT MET

11. Close the chamber door. The probe is now ready for use.

MET NOT MET

12. Record the HLD cycle on the log or printed sticker.

MET NOT MET

Notes:

- At the completion of a cycle, remove probe immediately to ensure faster warm up times. If probe remains in chamber, the Trophon will shut down heaters to ensure probe is protected. Therefore, warm up times may be longer.
- Sleep mode: to save power, Trophon will enter sleep mode after 2 hours of inactivity.
- Purging the Trophon of disinfectant is required if the device is to be moved or if Sonex-HL® has expired. Manufacturer's instructions for purging the Trophon must be meticulously followed.

NANOSONICS, INC

Trainer: _____

Trainee: _____

Date: 8-15-22

Appendix 6: TROPHON® COMPETENCY SHEET

NAME ALEXA TURK DATE 8/11/12

1. Put gloves on. Note: wear gloves at all times when required as described in the Trophon EPR user manual.

MET NOT MET

2. Pre-clean the probe before the high-level disinfection (HLD) cycle, following the probe manufacturer's instructions for use (IFUs).

MET NOT MET

3. Ensure the probe is clean and free of all visible debris, bioburden, gel, or other soil.

a. Use Sani Cloths to clean

b. Wipe the transducer cord and all surfaces of the transducer until it is visually clean. Use friction and work from cleanest to dirtiest areas.

c. Dry the transducer with a soft, dry cloth.

d. Visually inspect the transducer to ensure it is both clean and dry prior to HLD.

MET NOT MET

4. Load the clean, dry probe into the Trophon disinfection chamber ensuring:

a. The probe is secured high in the chamber with tip of probe above embossed line.

b. Probe does not contact the chamber wall at any point.

MET NOT MET

5. Place a new red Trophon chemical indicator (CI) into the indicator holder with the red side facing up.

a. Note: a new CI is to be used for every cycle.

MET NOT MET

6. Close the chamber door and confirm whether the probe is both clean and dry.

a. If yes, press Start.

b. If no, follow the LCD screen prompts.

MET NOT MET

7. At the end of the 7 minute HLD cycle, Trophon's LCD screen states: "CYCLE COMPLETE REMOVE AND WIPE PROBE."

MET NOT MET

8. Open the chamber door.

MET NOT MET

Cleaning, Disinfection, and Sterilization of Patient-Care Items

9. Remove CI, check CI color against the color chart on the CI carton and discard.
a. Note: BOTH the CI and LCD screen must indicate a successful cycle for the probe to be ready for use. If either the CI or Trophon LCD screen indicates a fail, repeat the cycle.

MET NOT MET

10. Remove and wipe the probe using a clean, dry single use cloth.

MET NOT MET

11. Close the chamber door. The probe is now ready for use.

MET NOT MET

12. Record the HLD cycle on the log or printed sticker.

MET NOT MET

Notes:

- At the completion of a cycle, remove probe immediately to ensure faster warm up times. If probe remains in chamber, the Trophon will shut down heaters to ensure probe is protected. Therefore, warm up times may be longer.
- Sleep mode: to save power, Trophon will enter sleep mode after 2 hours of inactivity.
- Purging the Trophon of disinfectant is required if the device is to be moved or if Sonex-HL® has expired. Manufacturer's instructions for purging the Trophon must be meticulously followed.

NANOSONICS, INC

Trainer: _____

Trainee: 8-15-22 _____

Date: 8-15-22 _____

Appendix 6: TROPHON® COMPETENCY SHEET

NAME SWANSON, Julie DATE 8/16/22

1. Put gloves on. Note: wear gloves at all times when required as described in the Trophon EPR user manual.

MET NOT MET

2. Pre-clean the probe before the high-level disinfection (HLD) cycle, following the probe manufacturer's instructions for use (IFUs).

MET NOT MET

3. Ensure the probe is clean and free of all visible debris, bioburden, gel, or other soil.

- a. Use Sani Cloths to clean
- b. Wipe the transducer cord and all surfaces of the transducer until it is visually clean. Use friction and work from cleanest to dirtiest areas.
- c. Dry the transducer with a soft, dry cloth.
- d. Visually inspect the transducer to ensure it is both clean and dry prior to HLD.

MET NOT MET

4. Load the clean, dry probe into the Trophon disinfection chamber ensuring:

- a. The probe is secured high in the chamber with tip of probe above embossed line.
- b. Probe does not contact the chamber wall at any point.

MET NOT MET

5. Place a new red Trophon chemical indicator (CI) into the indicator holder with the red side facing up.

- a. Note: a new CI is to be used for every cycle.

MET NOT MET

6. Close the chamber door and confirm whether the probe is both clean and dry.

- a. If yes, press Start.
- b. If no, follow the LCD screen prompts.

MET NOT MET

7. At the end of the 7 minute HLD cycle, Trophon's LCD screen states: "CYCLE COMPLETE REMOVE AND WIPE PROBE."

MET NOT MET

8. Open the chamber door.

MET NOT MET

Cleaning, Disinfection, and Sterilization of Patient-Care Items

9. Remove CI, check CI color against the color chart on the CI carton and discard.
a. Note: BOTH the CI and LCD screen must indicate a successful cycle for the probe to be ready for use. If either the CI or Trophon LCD screen indicates a fail, repeat the cycle.

MET NOT MET

10. Remove and wipe the probe using a clean, dry single use cloth.

MET NOT MET

11. Close the chamber door. The probe is now ready for use.

MET NOT MET

12. Record the HLD cycle on the log or printed sticker.

MET NOT MET

Notes:

- At the completion of a cycle, remove probe immediately to ensure faster warm up times. If probe remains in chamber, the Trophon will shut down heaters to ensure probe is protected. Therefore, warm up times may be longer.
- Sleep mode: to save power, Trophon will enter sleep mode after 2 hours of inactivity.
- Purging the Trophon of disinfectant is required if the device is to be moved or if Sonex-HL® has expired. Manufacturer's instructions for purging the Trophon must be meticulously followed.

NANOSONICS, INC

Trainer: _____

Trainee: _____

Date: 8-16-22

Appendix 6: TROPHON® COMPETENCY SHEET

NAME Leslie Aguirre DATE 8-15-11

- 1. Put gloves on. Note: wear gloves at all times when required as described in the Trophon EPR user manual.

MET NOT MET

- 2. Pre-clean the probe before the high-level disinfection (HLD) cycle, following the probe manufacturer's instructions for use (IFUs).

MET NOT MET

- 3. Ensure the probe is clean and free of all visible debris, bioburden, gel, or other soil.

- a. Use Sani Cloths to clean
- b. Wipe the transducer cord and all surfaces of the transducer until it is visually clean. Use friction and work from cleanest to dirtiest areas.
- c. Dry the transducer with a soft, dry cloth.
- d. Visually inspect the transducer to ensure it is both clean and dry prior to HLD.

MET NOT MET

- 4. Load the clean, dry probe into the Trophon disinfection chamber ensuring:

- a. The probe is secured high in the chamber with tip of probe above embossed line.
- b. Probe does not contact the chamber wall at any point.

MET NOT MET

- 5. Place a new red Trophon chemical indicator (CI) into the indicator holder with the red side facing up.

- a. Note: a new CI is to be used for every cycle.

MET NOT MET

- 6. Close the chamber door and confirm whether the probe is both clean and dry.

- a. If yes, press Start.
- b. If no, follow the LCD screen prompts.

MET NOT MET

- 7. At the end of the 7 minute HLD cycle, Trophon's LCD screen states: "CYCLE COMPLETE REMOVE AND WIPE PROBE."

MET NOT MET

- 8. Open the chamber door.

MET NOT MET

Cleaning, Disinfection, and Sterilization of Patient-Care Items

9. Remove CI, check CI color against the color chart on the CI carton and discard.
a. Note: BOTH the CI and LCD screen must indicate a successful cycle for the probe to be ready for use. If either the CI or Trophon LCD screen indicates a fail, repeat the cycle.

MET NOT MET _____

10. Remove and wipe the probe using a clean, dry single use cloth.

MET NOT MET _____

11. Close the chamber door. The probe is now ready for use.

MET NOT MET _____

12. Record the HLD cycle on the log or printed sticker.

MET NOT MET _____

Notes:

- At the completion of a cycle, remove probe immediately to ensure faster warm up times. If probe remains in chamber, the Trophon will shut down heaters to ensure probe is protected. Therefore, warm up times may be longer.
- Sleep mode: to save power, Trophon will enter sleep mode after 2 hours of inactivity.
- Purging the Trophon of disinfectant is required if the device is to be moved or if Sonex-HL® has expired. Manufacturer's instructions for purging the Trophon must be meticulously followed.

NANOSONICS, INC

Trainer: _____

Trainee: _____

Date: 8-15-21

Appendix 6: TROPHON® COMPETENCY SHEET

NAME D. Gunawan DATE 8-16-22

1. Put gloves on. Note: wear gloves at all times when required as described in the Trophon EPR user manual.

MET NOT MET

2. Pre-clean the probe before the high-level disinfection (HLD) cycle, following the probe manufacturer's instructions for use (IFUs).

MET NOT MET

3. Ensure the probe is clean and free of all visible debris, bioburden, gel, or other soil.

a. Use Sani Cloths to clean

~~b. Wipe the transducer cord and all surfaces of the transducer until it is visually clean. Use friction and work from cleanest to dirtiest areas.~~

c. Dry the transducer with a soft, dry cloth.

d. Visually inspect the transducer to ensure it is both clean and dry prior to HLD.

MET NOT MET

4. Load the clean, dry probe into the Trophon disinfection chamber ensuring:

a. The probe is secured high in the chamber with tip of probe above embossed line.

b. Probe does not contact the chamber wall at any point.

MET NOT MET

5. Place a new red Trophon chemical indicator (CI) into the indicator holder with the red side facing up.

a. Note: a new CI is to be used for every cycle.

MET NOT MET

6. Close the chamber door and confirm whether the probe is both clean and dry.

a. If yes, press Start.

b. If no, follow the LCD screen prompts.

MET NOT MET

7. At the end of the 7 minute HLD cycle, Trophon's LCD screen states: "CYCLE COMPLETE REMOVE AND WIPE PROBE."

MET NOT MET

8. Open the chamber door.

MET NOT MET

Cleaning, Disinfection, and Sterilization of Patient-Care Items

9. Remove CI, check CI color against the color chart on the CI carton and discard.
a. Note: BOTH the CI and LCD screen must indicate a successful cycle for the probe to be ready for use. If either the CI or Trophon LCD screen indicates a fail, repeat the cycle.

MET NOT MET

10. Remove and wipe the probe using a clean, dry single use cloth.

MET NOT MET

11. Close the chamber door. The probe is now ready for use.

MET NOT MET

12. Record the HLD cycle on the log or printed sticker.

MET NOT MET

Notes:

- At the completion of a cycle, remove probe immediately to ensure faster warm up times. If probe remains in chamber, the Trophon will shut down heaters to ensure probe is protected. Therefore, warm up times may be longer.
- Sleep mode: to save power, Trophon will enter sleep mode after 2 hours of inactivity.
- Purging the Trophon of disinfectant is required if the device is to be moved or if Sonex-HL® has expired. Manufacturer's instructions for purging the Trophon must be meticulously followed.

NANOSONICS, INC

Trainer: _____

Trainee: _____

Date: 8-16-22

Appendix 6: TROPHON® COMPETENCY SHEET

NAME UEDRANO, Daisy DATE 8/11/22

1. Put gloves on. Note: wear gloves at all times when required as described in the Trophon EPR user manual.

MET NOT MET

2. Pre-clean the probe before the high-level disinfection (HLD) cycle, following the probe manufacturer's instructions for use (IFUs).

MET NOT MET

3. Ensure the probe is clean and free of all visible debris, bioburden, gel, or other soil.

a. Use Sani Cloths to clean

b. Wipe the transducer cord and all surfaces of the transducer until it is visually clean. Use friction and work from cleanest to dirtiest areas.

c. Dry the transducer with a soft, dry cloth.

d. Visually inspect the transducer to ensure it is both clean and dry prior to HLD.

MET NOT MET

4. Load the clean, dry probe into the Trophon disinfection chamber ensuring:

a. The probe is secured high in the chamber with tip of probe above embossed line.

b. Probe does not contact the chamber wall at any point.

MET NOT MET

5. Place a new red Trophon chemical indicator (CI) into the indicator holder with the red side facing up.

a. Note: a new CI is to be used for every cycle.

MET NOT MET

6. Close the chamber door and confirm whether the probe is both clean and dry.

a. If yes, press Start.

b. If no, follow the LCD screen prompts.

MET NOT MET

7. At the end of the 7 minute HLD cycle, Trophon's LCD screen states: "CYCLE COMPLETE REMOVE AND WIPE PROBE."

MET NOT MET

8. Open the chamber door.

MET NOT MET

Cleaning, Disinfection, and Sterilization of Patient-Care Items

9. Remove CI, check CI color against the color chart on the CI carton and discard.
a. Note: BOTH the CI and LCD screen must indicate a successful cycle for the probe to be ready for use. If either the CI or Trophon LCD screen indicates a fail, repeat the cycle.

MET NOT MET

10. Remove and wipe the probe using a clean, dry single use cloth.

MET NOT MET

11. Close the chamber door. The probe is now ready for use.

MET NOT MET

12. Record the HLD cycle on the log or printed sticker.

MET NOT MET

Notes:

- At the completion of a cycle, remove probe immediately to ensure faster warm up times. If probe remains in chamber, the Trophon will shut down heaters to ensure probe is protected. Therefore, warm up times may be longer.
- Sleep mode: to save power, Trophon will enter sleep mode after 2 hours of inactivity.
- Purging the Trophon of disinfectant is required if the device is to be moved or if Sonex-HL® has expired. Manufacturer's instructions for purging the Trophon must be meticulously followed.

NANOSONICS, INC


Trainer: _____

Trainee: _____

Date: 8/15/22

trophon® EPR

TRAINING 2022

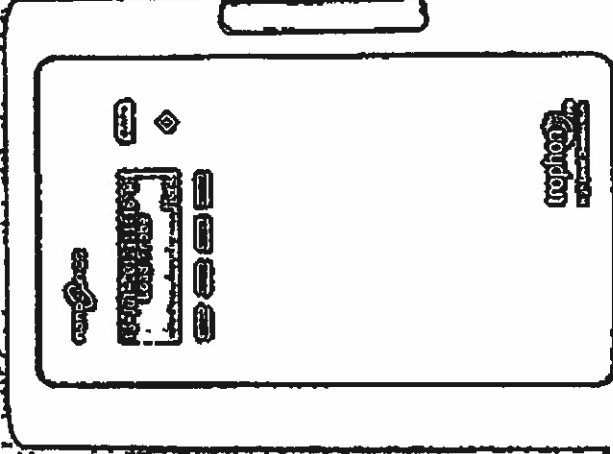

This certifies that

daisy medrano

has successfully completed the trophon EPR product training approved by Nanosonics Limited.



This certificate expires 12 months after the date of issue.

Date of issue: 12 Aug 2022

nanosonics
Infection Prevention. For Life.

trophon® EPR

TRAINING 2022

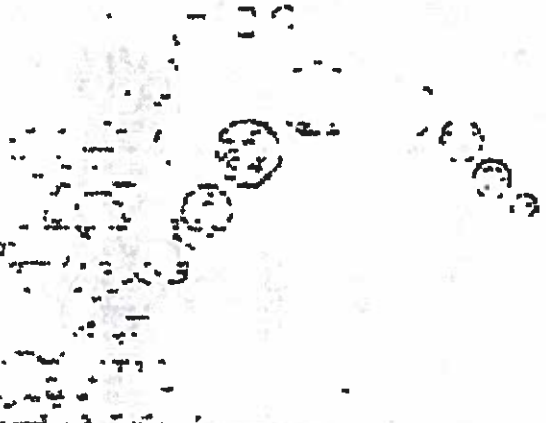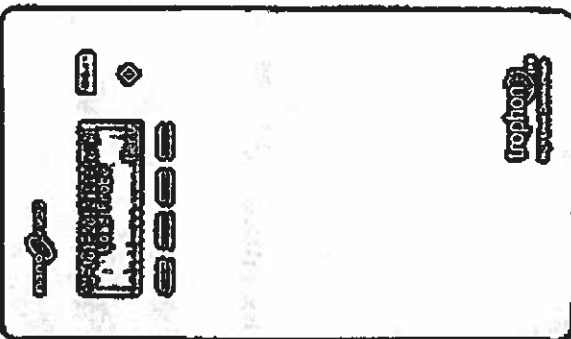

This certifies that

Andriy Martyniv

has successfully completed the trophon EPR product training approved by Nanosonics Limited.


This certificate expires 12 months after the date of issue.

Date of issue: 12 Aug 2022



trophon® EPR

TRAINING 2022

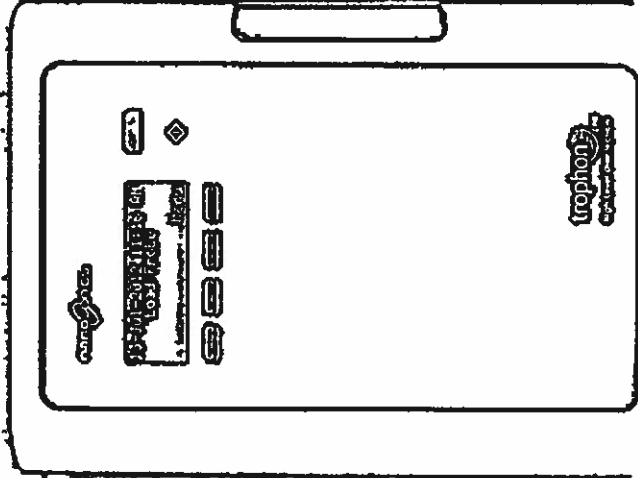

This certifies that

suhi Garcia

has successfully completed the trophon EPR product training approved by Nanosonics Limited.

This certificate expires 12 months after the date of issue.

Date of issue: 13 Aug 2022

nanosonics
Infection Prevention. For Life.

Notes:

At the completion of a cycle, remove probe immediately to ensure faster warm-up times.



If probe remains in chamber, the trophon EPR will shut down heaters to ensure probe is protected. Therefore, warm-up times may be longer.

Sleep mode – to save power, trophon EPR will enter sleep mode after two hours of inactivity.

Purging the trophon EPR of disinfectant is required if the device is to be moved or if NanoNebulant®/Sonex-HL® has expired. This involves removal of all disinfectant from the system. Purging may be implemented manually via trophon LCD screen (if relocating or transporting trophon) or it may be prompted automatically (if NanoNebulant®/Sonex-HL® cartridge has expired).

Empty the waste drawer when prompted by the trophon EPR LCD screen.

Note: Wear gloves when handling the waste drawer.

Trainer (Print Name):	<u>Megan Stacy, RPMS</u>
Trainer Signature:	
Title:	<u>Clin Apps Specialist</u>
Employee Signature:	
Date:	<u>8-15-22</u>

Nanosonics academy.com

trophon® EPR Competency

Name Gilbert Garcia

Date 8/15/22

1. TROPHON EPR OPERATION

1.1

Watch the trophon EPR training video (if not previously completed). This will take approximately 5 minutes.

Met Not Met



1.2

Complete and pass the online trophon EPR quiz.

Met Not Met



1.3

Check the trophon EPR Validated Probes List at your local Nanosonics/trophon website, under the trophon menu/trophon Probe compatibility. A hard copy of the list can be downloaded.

Met Not Met


4A	BL-D
4C-D	BL-RS
4C-RS	AB2-5
4C-9C	AB2-7
4D10L	AB2-7-D
4D18L	AB2-7-RS
4D3C-L	AD4-8
4D3C-R	AC2-5
4D8C	BE8C
4DE7C	BE9CS
4S	C1-5-D
4V-D	C1-5-RS
54R	C1-6-D
5C	C2-4-RC
5CS	C2-4-D
5S	C31
5S-02	C358
5S-4C	C36
5S-RS	C384
5S	C386
	C388
	C390
	C392
	C394
	C396
	C398
	C400
	C402
	C404
	C406
	C408
	C410
	C412
	C414
	C416
	C418
	C420
	C422
	C424
	C426
	C428
	C430
	C432
	C434
	C436
	C438
	C440
	C442
	C444
	C446
	C448
	C450
	C452
	C454
	C456
	C458
	C460
	C462
	C464
	C466
	C468
	C470
	C472
	C474
	C476
	C478
	C480
	C482
	C484
	C486
	C488
	C490
	C492
	C494
	C496
	C498
	C500
	C502
	C504
	C506
	C508
	C510
	C512
	C514
	C516
	C518
	C520
	C522
	C524
	C526
	C528
	C530
	C532
	C534
	C536
	C538
	C540
	C542
	C544
	C546
	C548
	C550
	C552
	C554
	C556
	C558
	C560
	C562
	C564
	C566
	C568
	C570
	C572
	C574
	C576
	C578
	C580
	C582
	C584
	C586
	C588
	C590
	C592
	C594
	C596
	C598
	C600
	C602
	C604
	C606
	C608
	C610
	C612
	C614
	C616
	C618
	C620
	C622
	C624
	C626
	C628
	C630
	C632
	C634
	C636
	C638
	C640
	C642
	C644
	C646
	C648
	C650
	C652
	C654
	C656
	C658
	C660
	C662
	C664
	C666
	C668
	C670
	C672
	C674
	C676
	C678
	C680
	C682
	C684
	C686
	C688
	C690
	C692
	C694
	C696
	C698
	C700
	C702
	C704
	C706
	C708
	C710
	C712
	C714
	C716
	C718
	C720
	C722
	C724
	C726
	C728
	C730
	C732
	C734
	C736
	C738
	C740
	C742
	C744
	C746
	C748
	C750
	C752
	C754
	C756
	C758
	C760
	C762
	C764
	C766
	C768
	C770
	C772
	C774
	C776
	C778
	C780
	C782
	C784
	C786
	C788
	C790
	C792
	C794
	C796
	C798
	C800
	C802
	C804
	C806
	C808
	C810
	C812
	C814
	C816
	C818
	C820
	C822
	C824
	C826
	C828
	C830
	C832
	C834
	C836
	C838
	C840
	C842
	C844
	C846
	C848
	C850
	C852
	C854
	C856
	C858
	C860
	C862
	C864
	C866
	C868
	C870
	C872
	C874
	C876
	C878
	C880
	C882
	C884
	C886
	C888
	C890
	C892
	C894
	C896
	C898
	C900
	C902
	C904
	C906
	C908
	C910
	C912
	C914
	C916
	C918
	C920
	C922
	C924
	C926
	C928
	C930
	C932
	C934
	C936
	C938
	C940
	C942
	C944
	C946
	C948
	C950
	C952
	C954
	C956
	C958
	C960
	C962
	C964
	C966
	C968
	C970
	C972
	C974
	C976
	C978
	C980
	C982
	C984
	C986
	C988
	C990
	C992
	C994
	C996
	C998
	C1000

1.4

Put gloves on.

Note: Wear gloves at all times when required as described in the *trophon® EPR User Manual*.


Met Not Met

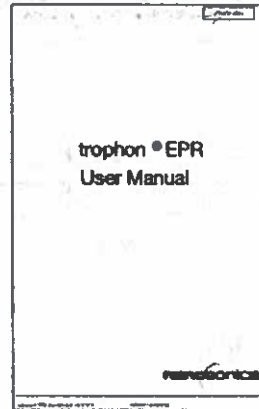

1.5

Pre-clean the probe before the high level disinfection (HLD) cycle, following the probe manufacturer's instructions.

Met Not Met

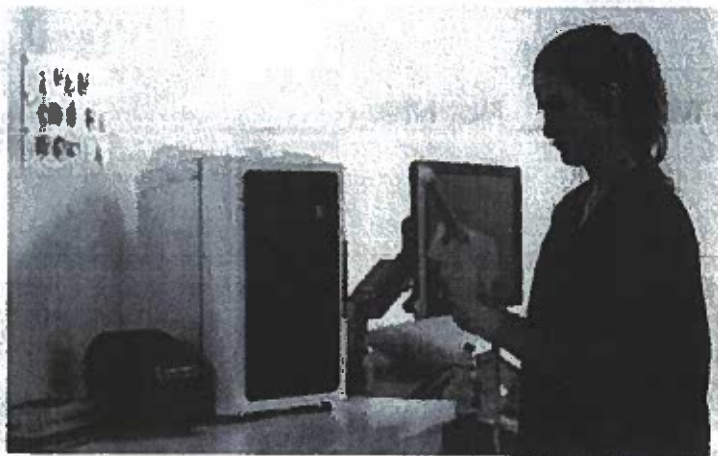


1.6

Ensure the probe is clean and free of all visible debris, bioburden, gel or other soil.

Met Not Met

Cleaning Process

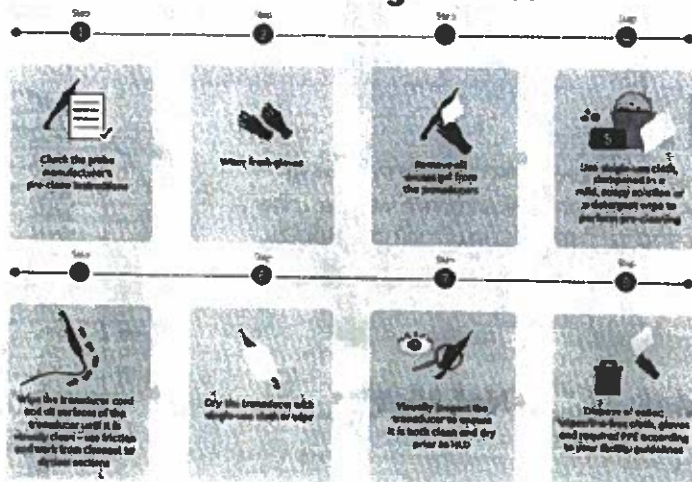

1.7

Dry the probe – ensure there is no visible moisture or drops of liquid on the probe.

Met Not Met



1.8

Load the clean, dry probe into the trophon EPR disinfection chamber, ensuring:

- Probe is secured high in the chamber with tip of probe above embossed line
- Probe does not contact the chamber wall at any point

If you have a curved probe, check the *trophon EPR Validated Probes List* at your local Nanosonics/trophon website or download the list from the Nanosonics training site, to see if you need the CPP accessory.

Met Not Met



1.9

Place a Chemical Indicator into the locator at the base of the chamber door, red side up.

Note: A new Chemical Indicator to be used for every cycle

Met Not Met



1.10

Close the chamber door and confirm whether the probe is both clean and dry

- If yes, select "Yes" then press Start
- If no, follow the LCD screen prompts

Met Not Met



1.11

At the end of the seven minute HLD cycle, trophon EPR's LCD screen states: CYCLE COMPLETE REMOVE AND WIPE PROBE.

Met Not Met



1.12


Wear new gloves

Met Not Met

1.13

Open the chamber door.

Met Not Met

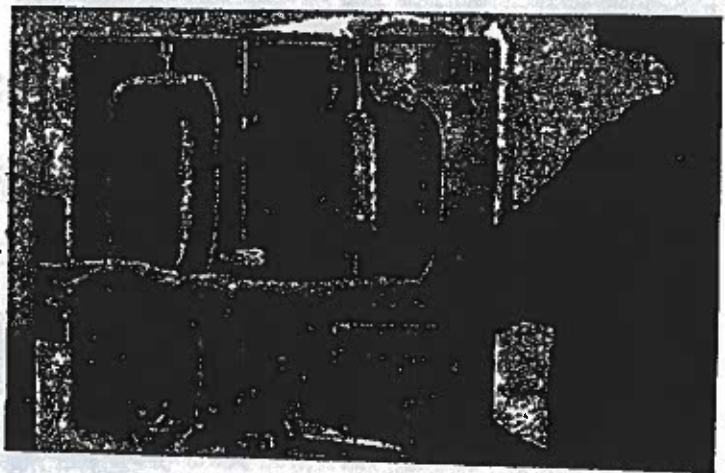


1.14

Remove CI, check CI colour against the colour chart on the CI box and discard.

Note: BOTH the CI and LCD screen must indicate a successful cycle for the probe to be ready for use. If either the CI or Trophon LCD screen indicates a fail, repeat the cycle.

Met Not Met

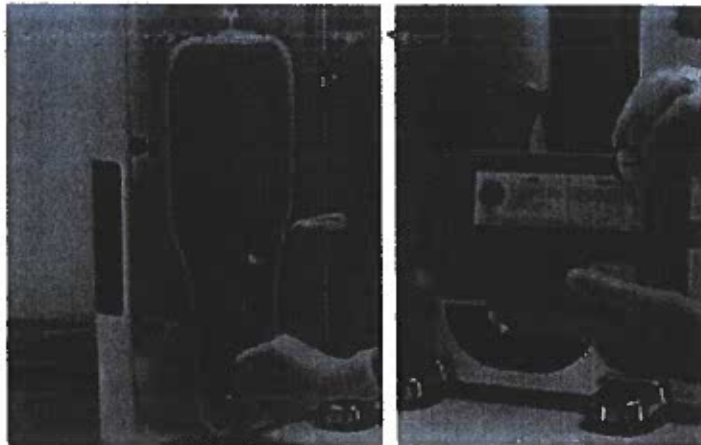

1.15

Remove and wipe the probe using a dry, single use cloth.

Met Not Met



1.16

Insert the probe into a trophon Clean Ultrasound Probe Cover and wrap the securing tie around the cable strain relief.



Met Not Met



1.17

Record the high level disinfection cycle on the log or printed sticker. Store the probe correctly for later use

Met Not Met

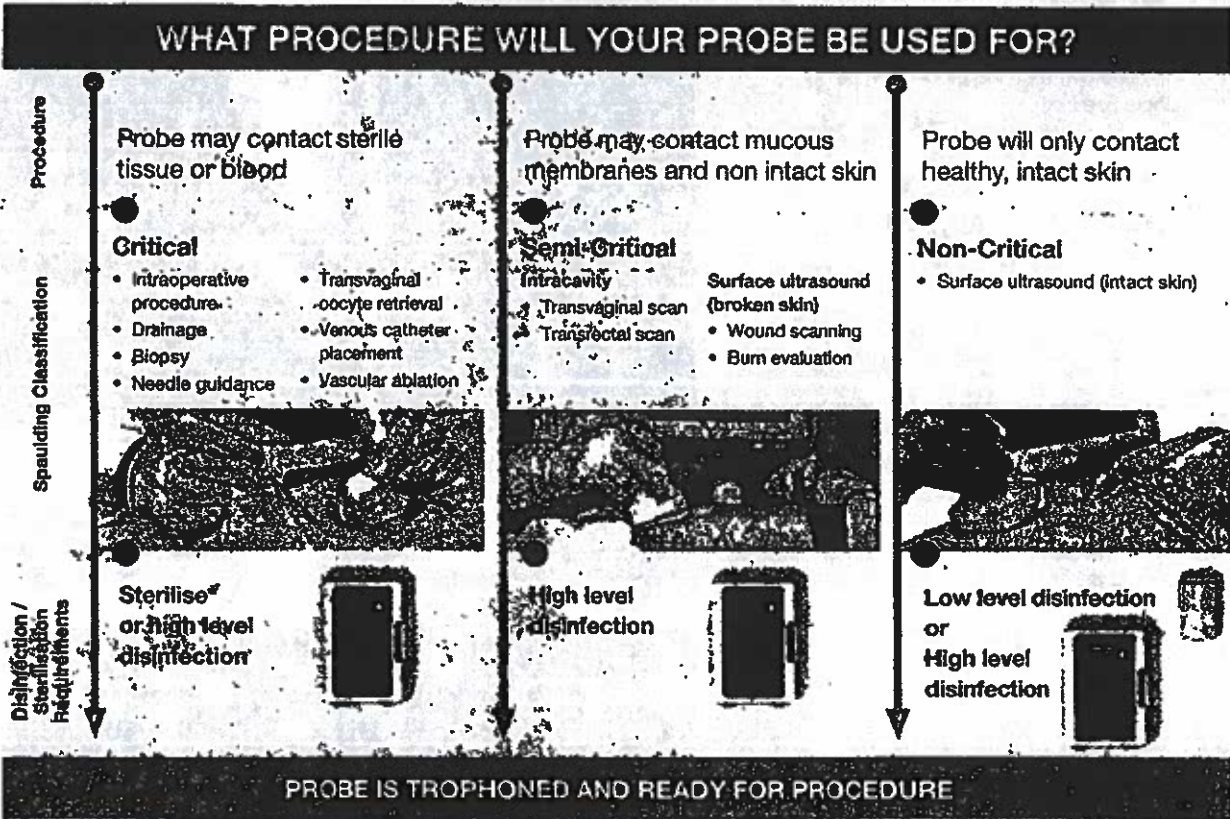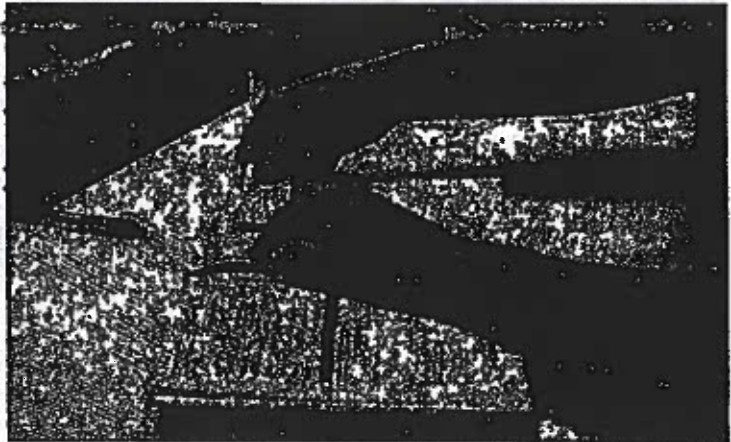

* Critical probes should be sterilised, or can also be high level disinfected and used with a sterile sheath.¹

¹ Rutala WA, Weber DJ, HICPAC. Guideline for Disinfection and Sterilization in Healthcare Facilities. USA: Centers for Disease Control; CDC 2008 (<https://www.cdc.gov/infectioncontrol/pd/guidelines/disinfectionguidelines.pdf>)

trophon® EPR Competency

Name Arnold Sabater RN

Date 8/15/22

1. TROPHON EPR OPERATION

1.1

Watch the trophon EPR training video (if not previously completed). This will take approximately 5 minutes.

Met Not Met



1.2

Complete and pass the online trophon EPR quiz.

Met Not Met




1.3

Check the trophon EPR Validated Probes List at your local Nanosonics/trophon website, under the trophon menu/trophon Probe compatibility. A hard copy of the list can be downloaded.

Met Not Met

40-A	8L-D
4C-D	8L-RS
4C-RS	AB2-6
4C-SC	AB2-7
4D10L	AB2-7-D
4D18L	AB2-7-RS
4D2C-L	AB4-8
4D6C-L	AC2-5
4D8C	BE9C
4DE7C	BE9CS
4S	C1-4-D
4V-D	C1-5-RS
4W lens	C1-6-D
5C	C2-6-RC
5C lens	C2-9-D
5S	C31
5S-62	C368
5S-LC	C36
5S-RS	C364
6S	C366
	C4
	8L-D
	8L-RS
	AB2-6
	AB2-7
	AB2-7-D
	AB2-7-RS
	G35
	712L
	112L-RS
	1738
	1738-RS
	ICS-9
	ICS-9-D
	ICS-9H
	CS-8W-RS
	L3-12-D
	L39
	L4-12

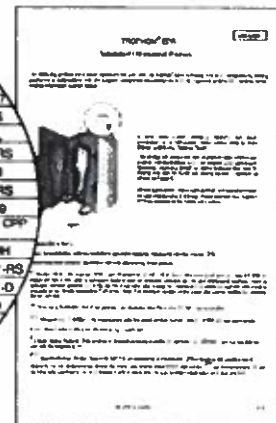


1.4

Put gloves on

Note: Wear gloves at all times when required as described in the *trophon® EPR User Manual*.

Met Not Met

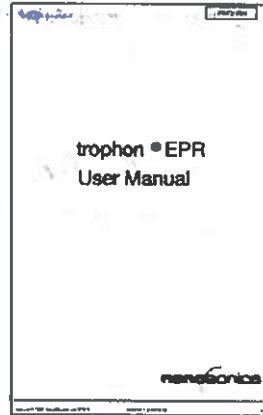

1.5

Pre-clean the probe before the high level disinfection (HLD) cycle, following the probe manufacturer's instructions.

Met Not Met




1.6

Ensure the probe is clean and free of all visible debris, bioburden, gel or other soil.

Met Not Met

Cleaning Process

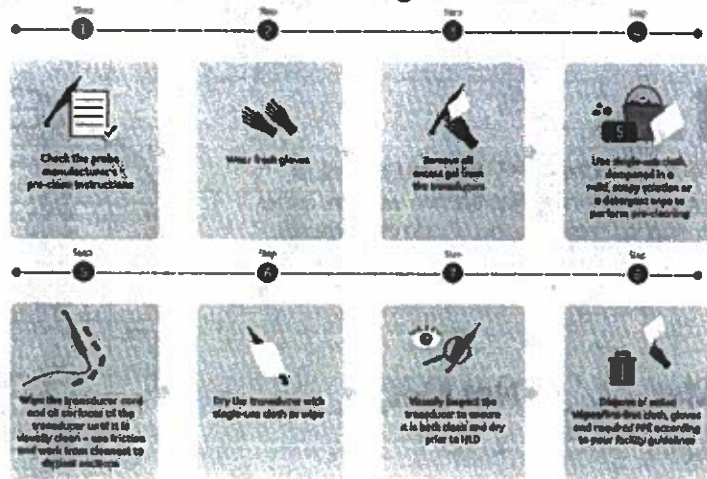

1.7

Dry the probe – ensure there is no visible moisture or drops of liquid on the probe.

Met Not Met



1.8

Load the clean, dry probe into the trophon EPR disinfection chamber, ensuring:

- Probe is secured high in the chamber with tip of probe above embossed line
- Probe does not contact the chamber wall at any point

If you have a curved probe, check the *trophon EPR Validated Probes List* at your local Nanosonics/trophon website or download the list from the Nanosonics training site, to see if you need the CPP accessory.

Met Not Met



1.9

Place a Chemical Indicator into the locator at the base of the chamber door, red side up.

Note: A new Chemical Indicator to be used for every cycle

Met Not Met



1.10

Close the chamber door and confirm whether the probe is both clean and dry

- If yes, select "Yes" then press Start
- If no, follow the LCD screen prompts

Met Not Met



1.11

At the end of the seven minute HLD cycle, trophon EPR's LCD screen states: CYCLE COMPLETE REMOVE AND WIPE PROBE.

Met Not Met



1.12

Wear new gloves.

Met Not Met

1.13

Open the chamber door.

Met Not Met




1.14

Remove CI, check CI colour against the colour chart on the CI box and discard.

Note: BOTH the CI and LCD screen must indicate a successful cycle for the probe to be ready for use. If either the CI or Trophon LCD screen indicates a fail, repeat the cycle.

Met Not Met

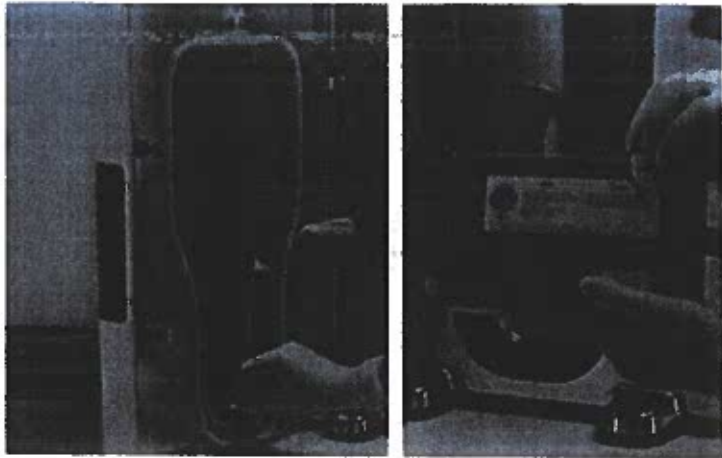

1.15

Remove and wipe the probe using a dry, single use cloth.

Met Not Met



1.16

Insert the probe into a trophon Clean Ultrasound Probe Cover and wrap the securing tie around the cable stain relief.

Met Not Met



7
 e high level disinfection cycle on the
 ted sticker. Store the probe correctly
 se.



Not Met



WHAT PROCEDURE WILL YOUR PROBE BE USED FOR?

Probe may contact sterile tissue or blood

Critical

- Intraoperative procedure
- Drainage
- Biopsy
- Needle guidance
- Transvaginal oocyte retrieval
- Venous catheter placement
- Vascular ablation



**Sterilise⁺
or high level
disinfection**



Probe may contact mucous membranes and non intact skin

Semi-Critical

- Intracavity**
- Transvaginal scan
- Transrectal scan
- Surface ultrasound (broken skin)**
- Wound scanning
- Burn evaluation



**High level
disinfection**



Probe will only contact healthy, intact skin

Non-Critical

- Surface ultrasound (intact skin)



**Low level disinfection
or
High level
disinfection**



PROBE IS TROPHONED AND READY FOR PROCEDURE

Probes should be sterilised, or can also be high level disinfected and used with a sterile sheath.¹

¹A, Weber DJ, HICPAC. Guideline for Disinfection and Sterilization in Healthcare Facilities. USA: Centers for Disease Control; B (<https://www.cdc.gov/infectioncontrol/pdf/guidelines/disinfectionguidelines.pdf>)

Notes:

At the completion of a cycle, remove probe immediately to ensure faster warm-up times.

If probe remains in chamber, the trophon EPR will shut down heaters to ensure probe is protected. Therefore, warm-up times may be longer.

Sleep mode – to save power, trophon EPR will enter sleep mode after two hours of inactivity.

Purging the trophon EPR of disinfectant is required if the device is to be moved or if NanoNebulant®/Sonex-HL® has expired. This involves removal of all disinfectant from the system. Purging may be implemented manually via trophon LCD screen (if relocating or transporting trophon) or it may be prompted automatically (if NanoNebulant®/Sonex-HL® cartridge has expired).

Empty the waste drawer when prompted by the trophon EPR LCD screen.

Note: Wear gloves when handling the waste drawer.

Trainer (Print Name):	<i>Megan Stacy, RDMS</i>
Trainer Signature:	[Redacted]
Title:	<i>Clin Apps Specialist</i>
Employee Signature:	[Redacted]
Date:	<i>8/15/22</i>

Appendix 5: TROPHON® Log

High-Level Disinfection Log
 For Trophon EPR Serial # 18640-009

START A NEW PAGE WHEN LOADING A NEW SONEX CARTRIDGE OR NEW BOX OF CHEMICAL INDICATORS

CARTRIDGE REPLACEMENT: WHEN LOADING A NEW SONEX CARTRIDGE, START A NEW PAGE. DOCUMENT DATE LOADED, LOT # & EXP. DATE.

CHEMICAL INDICATOR BOX: WHEN OPENING A NEW BOX OF CHEMICAL INDICATORS, START A NEW PAGE. DOCUMENT LOT # AND EXP. DATE.

CYCLE RESULTS: PLACE LABEL IN SQUARES BELOW.

LABELING YOUR READY-TO-USE PROBE: INITIAL IN THE SMALL BOX BELOW ("INIT") AND LABEL IS PLACED ON THE CLEAN US PROBE COVER.

CARTRIDGE REPLACEMENT		CHEMICAL INDICATOR REPLACEMENT	
Date Loaded & Lot #	Expiration Date	Lot Batch #	Expiration Date
8/18/22 C11001	2023-10-06	21 P08 10	2023-03-02

trophon
 17/08/2022 07:58 SN:18640-009
 Disinfection: PASS Cycle #: 5490
 Indicator: PASS FAIL
 Operator: [Signature]
 Probe: Red
 Notes: Test

trophon
 17/08/2022 08:48 SN:18640-009
 Disinfection: PASS Cycle #: 5491
 Indicator: PASS FAIL
 Operator: [Signature]
 Probe: Green
 Notes: Test Run

2 nd Label Placed on Clean Ultrasound Probe Cover	2 nd Label Placed on Clean Ultrasound Probe Cover	2 nd Label Placed on Clean Ultrasound Probe Cover
2 nd Label Placed on Clean Ultrasound Probe Cover	2 nd Label Placed on Clean Ultrasound Probe Cover	2 nd Label Placed on Clean Ultrasound Probe Cover
2 nd Label Placed on Clean Ultrasound Probe Cover	2 nd Label Placed on Clean Ultrasound Probe Cover	2 nd Label Placed on Clean Ultrasound Probe Cover
2 nd Label Placed on Clean Ultrasound Probe Cover	2 nd Label Placed on Clean Ultrasound Probe Cover	2 nd Label Placed on Clean Ultrasound Probe Cover

American Health Care Centers

List of Staff Competent to Use Trophon EPR

Name	Initial
Arnold Sabater RN (super trainer) -----	[Redacted]
Andrei MARTyniv, SA (super trainer) -----	[Redacted]
Gilbert Garcia, ORT (super trainer) -----	GG
Daisy Medrano, ORT -----	D.M.
Julie Swanson, Faculty Administrator -----	JSS
Devi Gunawan, APN -----	DG
Leslie Aguirre, ULS tech -----	LA
Alexa Tulea, OR staff -----	AT
Marie Frukacz, Manager -----	MF
Natasha Kukuriza, ULS tech -----	_____ pending competency

[Redacted] Red = SN: 64400530

[Redacted] Green = SN: B7626609